VERNEBLER NEBUFIRST® / NEBUJUNIOR®

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DE

PIKTOGRAMM-LEGENDE

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Bestellnummer
















Einpatientenprodukt

Inhaltsangabe in Stück

Hersteller

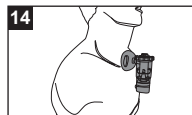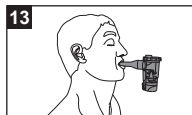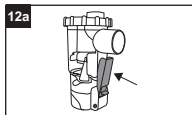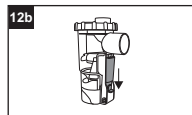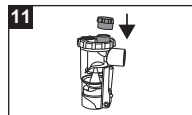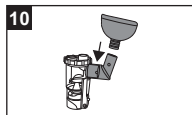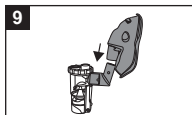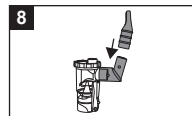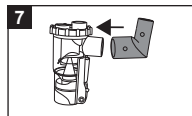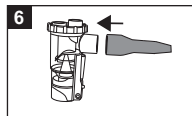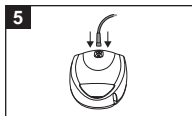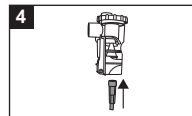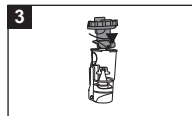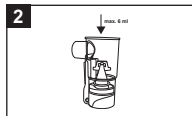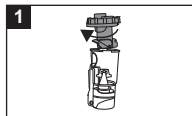I. VORWORT

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NEBULISER NEBUFIRST® / NEBUJUNIOR®

I. FOREWORD

These brief instructions for use are valid for the NEBUFIRST® / NEBUJUNIOR® nebuliser. The instructions for use serve to inform the physician, nursing staff and patient/user in order to ensure technically correct handling.

Please read the instructions for use carefully before using the product for the first time!

Please also observe the instructions for use of the NEBUFIRST® / NEBUJUNIOR® INHALATOR.

II. INTENDED USE

Humidification of breathing air through inhalation of steam and medications and other inhalants approved for use in inhalation therapy.

III. IMPORTANT SAFETY INSTRUCTIONS

CAUTION

Do not perform inhalation in the airways until you have been briefed by specialist staff.

Only use the device and accessories if they are clean and disinfected.

Unless otherwise prescribed by the physician, inhalation should be performed 2 – 3 times daily for approx. 10 – 15 minutes.

Once phlegm has loosened, it should be coughed up, or in the case of tracheostomy patients, suctioned using a suctioning device and a sterile suction catheter.

Draughts and cold temperatures should be avoided immediately after inhalation.

IV. INSTRUCTIONS

1. CONNECTING THE NEBULISER

- Pull off the top section of the nebuliser by twisting both parts (anti-clockwise) gently against each other ①
- In the lower section of the nebuliser ②, pour in the exact quantity of the medication/inhalant in accordance with the instructions provided by the manufacturer or the attending physician.
- Observe the maximum capacity of 6 ml.
- Place the top section of the nebuliser back onto the lower section. Connect the nebuliser by twisting both parts (clockwise) gently against each other. Ensure that these fit together securely ③.

2. CONNECTING THE AIR TUBE

- Check whether the air tube is bent or kinked. Connect the air tube ④ with the lower opening of the nebuliser.
- Then connect the air tube with the air connection of the device ⑤.

3. CONNECTING ATTACHMENTS

- Attach the mouthpiece ⑥ to the side connection of the nebuliser. Please also use the adapter ⑦ included in the delivery in order to make use of the nosepiece ⑧, the inhalation mask ⑨ or the LARYVOX® tracheostoma mask ⑩.
- In order to increase breathing resistance, you can place the PEP cap included in the delivery on to the expiration valve on the top section of the nebuliser ⑪.

4. NEBULISER

- Hold the nebuliser in an upright position where possible.
- Two applications are possible: continuous inhalation and temporary inhalation. For continuous inhalation, push the locking device (12 b) downwards in the direction indicated by the arrow in order to secure the blue lever for the desired duration of nebulisation. For temporary inhalation, hold down the blue lever (12 a).
- To perform inhalation, sit in a relaxed pose and breath in and out evenly.
- To allow your airways to become accustomed to the aerosol, do not breathe in too deeply during the first breaths you inhale.
- In order to increase breathing resistance, you can place the PEP cap on to the expiration valve on the top section of the nebuliser (13).

5. ATTACHMENTS

- When using the mouthpiece, enclose it fully with your lips (14) the mouthpiece must be connected via the adapter.
- During inhalation, the LARYVOX® tracheostoma mask (special accessories) must be connected to the nebuliser using the adapter.
- Press the inhalation mask/tracheostoma mask gently over your face/tracheostomy (15), relax, and breath in.
- When breathing out, the inhalation mask/tracheostoma mask should be held a small distance (the width of a finger) away from your face/tracheostomy.

V. CLEANING, DISINFECTION, REPROCESSING

Please refer to the instructions for use of the NEBUFIRST® / NEBUJUNIOR® INHALATOR for guidelines on cleaning, disinfection and reprocessing.

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Medizintechnik-Vertrieb GmbH