



- ▶ Speaking valve - transparent, lightweight
- ▶ With silicone membrane
- ▶ Membrane closed in basic position (“closed-position valve”) – opens during inspiration.
 Advantage: No enhanced exhalation pressure required to close the valve and to facilitate speaking. Less fatigue for the patient, better speaking quality
- ▶ With multiadapter: fits to both 15 mm connector (UNI) and 22 mm adapter (KOMBI)
- ▶ To be used by tracheotomized patients, not to be used by laryngectomized patients
- ▶ To be used only and exclusively in combination with a phonation cannula.
- ▶ COMBIPHON® O₂ with oxygen supply port
- ▶ Sterile
- ▶ Individually packed

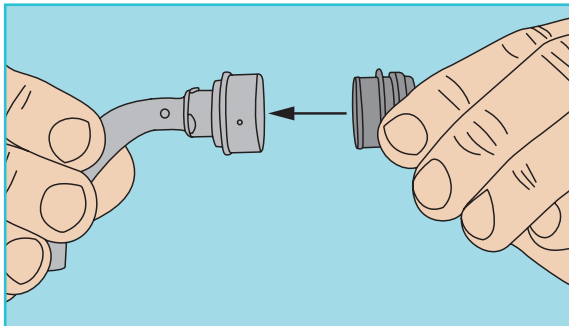

STERILE EO

Sterilized using ethylene oxid

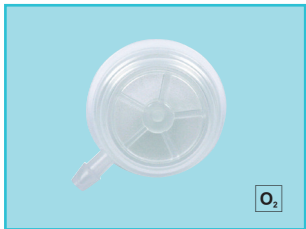

Do not resterilize



Keep away from sunlight



Keep dry



Single patient use



Consult instructions for use



CE marking identification number of the notified body

Product made in Germany

Order information:

COMBIPHON®

COMBIPHON® O₂

REF 27131

REF 27132