


DURACUFF® VARIO-KOMBI LINGO-PHON SUCTION 1 ICV + 1 ICFK + 1 COMBIPHON®

- ▶ Tracheostomy tube with two inner cannulas: one with a rotating 15 mm swivel connector (ICV / VARIO) to buffer external forces like e.g tension and torque forces exerted by the ventilation tubing – one sieve-fenestrated with a 22 mm adapter (ICFK / KOMBI) to connect a speaking valve
- ▶ With low-pressure cuff (CUFF) and integrated suction line (SUCTION)
- ▶ With suction line (SUCTION), integrated into the outer cannula to retain a smooth and atraumatic surface of the tracheostomy tube to avoid potential irritations around the tracheostoma
- ▶ With DURACUFF® Clip for locking the inner cannula into the outer
- ▶ With attachable COMBIPHON® closed-position speaking valve
- ▶ Soft, pliable, lightweight
- ▶ Transparent, inconspicuous design
- ▶ Set includes broad and adjustable tube holder (OPTIFLAUSCH®), obturator and product ID card
- ▶ Three lengths: Standard length, medium (M) and extra-long (XL)

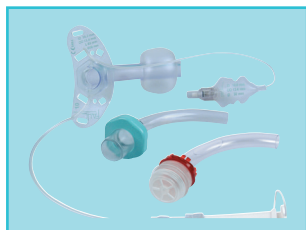

Order information:

DURACUFF® VARIO-KOMBI LINGO-PHON SUCTION,
1 ICV + 1 ICFK + 1 COMBIPHON®, SIZE 7-12



REF 13042

DURACUFF® VARIO-KOMBI LINGO-PHON SUCTION M,
1 ICV + 1 ICFK + 1 COMBIPHON®, SIZE 7-12



REF 13142

DURACUFF® VARIO-KOMBI LINGO-PHON SUCTION XL,
1 ICV + 1 ICFK + 1 COMBIPHON®, SIZE 7-12



REF 13242

When ordering, please state the size after the item number!

- STERILE EO** Sterilized using ethylene oxide
- Do not resterilize
- Keep away from sunlight
- Single patient use
- Consult instructions for use

Rx only CAUTION: United States Federal law restricts this device to sale, distribution and use by or on order of a physician or a licensed practitioner (in the US only).

CE 0482 CE marking identification number of the notified body

Product made in Germany

Size	O.D. neck flange mm	O.D. Tip mm	I.D. mm O-cannula	I.D. mm I-cannula	Length mm standard	Length mm medium (M)	Length mm extra long (XL)
7	10.2	8.8	7.0	5.5	65	75	100
8	11.2	9.8	8.0	6.5	70	80	105
9	12.4	11.2	9.0	7.5	75	85	105
10	13.4	12.2	10.0	8.5	85	95	110
11	14.4	13.2	11.0	9.5	87	97	115
12	15.4	14.2	12.0	10.2	90	100	120

O.D. Tip = Outer diameter of the outer tube at the tip of the cannula · O.D. neck flange = Outer diameter behind the neck flange · I.D. = Inner diameter at the tip of the cannula · Length = Length over the outer curve