







FAHL® TRACHEO-SAFE

- ▶ Tracheostoma place holder to be inserted into the tracheostoma and to secure the airway while minimizing the contact surface with the delicate tracheal wall
- ▶ Made of soft silicone
- ▶ Consists of a cannula that, at its lower end, runs into an elongated half-tube, adjusted to match the trachea and fitting anatomically against the wall of the trachea
- ▶ Neck flange with two lateral openings to attach a tube holder
- ▶ Adjustable neck flange to individually determine the space between the neck surface and the trachea
- ▶ A detachable sealing plug is located at the upper (distal) end of the place holder which sticks out of the trachea to seal the tracheostoma when weaning a patient off the tracheostomy tube, e.g. after long-term ventilation
- ▶ By removing the plug it is possible to suction tracheal secretions through the cannula
- ▶ Available in two different lengths, 45 and 70 mm

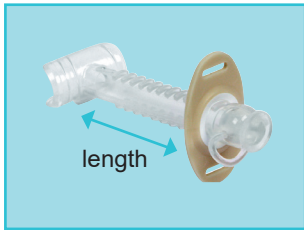

 Sterilized using ethylene oxid
 Do not re-sterilize
 Keep away from sunlight
 Single patient use
 Consult instructions for use
CE 0482 CE marking identification number of the notified body
Product made in Germany

Order information:

FAHL® Tracheo-Safe C, 45 mm length

REF 20895

FAHL® Tracheo-Safe P, 70 mm length

REF 20898