


# LARYVOX® EMERGENCY-SET LARYNGECTOMY PRO

Complete Set for the acute care of laryngectomized patients in case of emergency

- ▶ 1x LARYVOX® SOS safe-mask  
Emergency ventilation mask for laryngectomees
- ▶ 1x TRACHEOTEC® VARIO CUFF  
Tracheostomy tube with cuff, size 07
- ▶ 1x LARYVOX® Tracheospreizer (tracheal dilator)  
Using the LARYVOX® Tracheospreizer, the cannula can be inserted safely and as quickly as possible
- ▶ 1x LARYVOX® Connect  
Connects the 22 mm adapter (KOMBI) of the SOS safe-mask with a 15 mm connector (UNI) of e.g. a tracheostomy tube
- ▶ 5x LARYVOX® O<sub>2</sub> Connect  
With multiadapter; to be connected to the Laryvox Tape base plate with 22 mm adapter (KOMBI) or to a tracheostomy tube with 15 mm (UNI) connector. The oxygen supply port allows the administration of extra oxygen
- ▶ 5x LARYVOX® Tape COMFORT / XL Oval  
Base plate with strong adhesive strength and an extra-large surface, with 22 mm adapter (KOMBI)
- ▶ 3x OPTIFLUID®  
Stoma oil towels with colourless neutral oil to lubricate the Tracheotec tracheostomy tube before insertion
- ▶ 1x Oxygen connection tube
- ▶ 1x Single-use syringe / 10ml

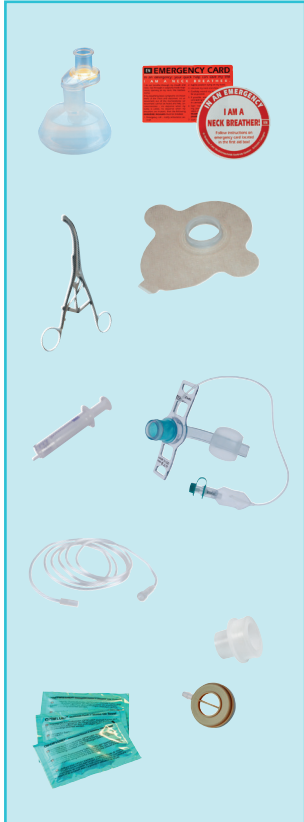

**For all components of this kit separate and detailed data sheets are available.**

## Order information:

**LARYVOX® EMERGENCY-SET Laryngectomy PRO**

**REF 48007**