

LARYVOX® ADHESIVE FOAM DISC


CUSHION RING

- ▶ LARYVOX® Adhesive Foam Disc is a foam ring with a central circular opening
- ▶ The front and back are provided with an adhesive surface covered by a protective film
- ▶ It is used to align the skin level of an irregular tracheostoma that is situated below the skin level and thus creates a flatter surface for the fixation of a base plate
- ▶ In combination with LARYVOX® Tape or a sealing plate, the contact surface of the tape is optimised, thus enabling a longer wearing time
- ▶ The speaking option is improved by exclusion of undesired losses of air
- ▶ The different size versions small (REF 48460) and large (REF 48461) serve for individual adaptation in order to achieve an even skin level to be covered with LARYVOX® Tape



REF 48460

CE

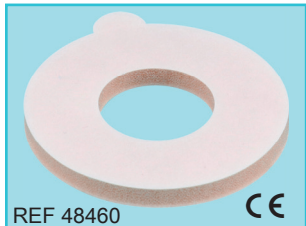


REF 48460

CE



REF 48461

CE

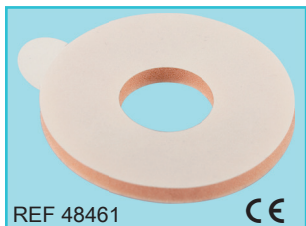

REF 48461

CE

Order information:

LARYVOX® ADHESIVE FOAM DISC small, Pack of 30

REF 48460

LARYVOX® ADHESIVE FOAM DISC large, Pack of 30

REF 48461