

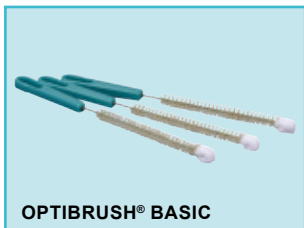


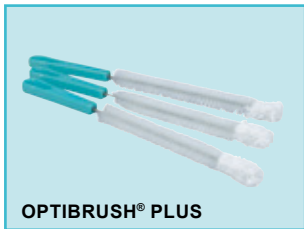


- ▶ Special brush for gentle cleaning of tracheostomy tube
- ▶ Flexible
- ▶ Bristles in dense array around the stainless steel wire base which is twisted to avoid sharp edges which might damage the tracheostomy tube
- ▶ OPTIBRUSH® (REF 31850) comes with a rounded tip to avoid scratching of the soft tracheostomy tube material
- ▶ The brush tip of OPTIBRUSH® PLUS (REF 31855) is surrounded by a soft, wipe-active and absorbent fibre top which further reduces the risk to damage the tube and which ensures a more effective cleaning
- ▶ The brush tip of OPTIBRUSH® Basic (REF 31800) is surrounded by a soft, approx. 1 cm wide, wipe-active and absorbent wool tuft which further reduces the risk to damage the tube and which ensures a more effective cleaning
- ▶ Wide size range

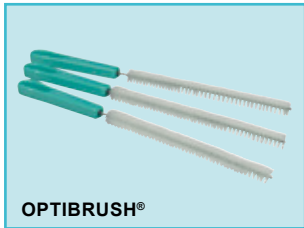

 Single patient use
 Consult instructions for use
 CE marking

Order information:	
OPTIBRUSH® , Size 4 - 14, Pack of 4	REF 31850
OPTIBRUSH® PLUS , Size 5 - 14, Pack of 4	REF 31855
OPTIBRUSH® Basic , Size 6 - 14, Pack of 4	REF 31800