


TRACHLINE® UNI LINGO-PHON CUFF 1 ICS + 1 ICU + 1 ICSU + SV


- ▶ Tracheostomy tube with CUFF for long-term use, max. service life up to 29 days
- ▶ Made of thermosensitive medical synthetic materials
- ▶ Soft, pliable, lightweight
- ▶ Transparent, inconspicuous design
- ▶ Rounded tube tip to prevent irritation of the mucous membranes
- ▶ Bayonet locks enable simple, secure attachment of the inner cannula and accessories
- ▶ Set consists of a sieved inner cannula (ICS), a sieved inner cannula with 15 mm connector (ICSU), an inner cannula with 15 mm connector (ICU), a 15 mm connector (UNI), a cough cap, a sealing cap, a locking ring, a speaking valve, an adjustable tube holder, an obturator as well as a cannula passport
- ▶ TRACHLINE® PHON VALVE (speaking valve) and FAHL® UNI SEAL (sealing cap) are also available separately



CE 0482

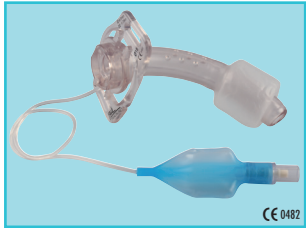

CE 0482

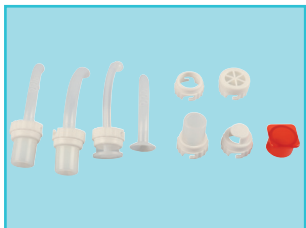

Order information:

TRACHLINE® UNI LINGO-PHON CUFF Size 5-10
1 ICS + 1 ICU + 1 ICSU + SV



REF 13801

FAHL® UNI SEAL

REF 13830

TRACHLINE® PHON VALVE

REF 27141

When ordering, please state the size after the item number!

Size	O.D. neck flange mm	I.D. O.C. mm	I.D. I.C. mm	Angel ca.	Cuff (Ø) mm	Length mm
5	10,6	6.5	5.0	115°	18	62
6	11,7	7.5	6.0	110°	22	70
7	12,8	8.5	7.0	105°	22	72
8	13,8	9.5	8.0	100°	26	80
9	14,8	10.5	9.0	95°	30	82
10	15,8	11.5	10.0	90°	32	84

O.D. neck flange = Outer diameter behind the neck flange · I.D. = Inner diameter at the tip of the cannula · O.C. = Outer cannula · I.C. = Inner cannula · Length = Length over middle curve