
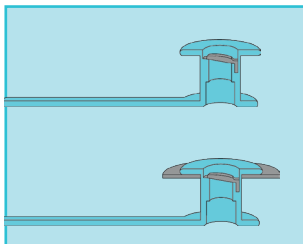
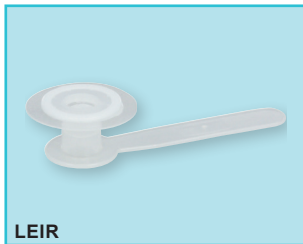


# BLOM-SINGER® CLASSIC™ INDWELLING IR / LEIR SPECIAL ORDER (INCREASED RESISTANCE / LARGE ESOPHAGEAL FLANGE) VOICE PROSTHESIS

- ▶ The Blom-Singer® Special Order Series offers a large assortment of special indwelling prostheses that cover specific requirements
- ▶ All variants are made of medical grade silicone and are furnished with specific product features to suit the respective indications
- ▶ Blom-Singer® CLASSIC™ Indwelling IR (Increased Resistance) is furnished with a special high-resistance valve
- ▶ This valve is equipped with a reinforced silicone flap hinge and prevents unintentional opening of the valve in the event of excessively high tracheo-oesophageal suction pressure
- ▶ Optionally, this version is also available as a „Large Flange“ version: Blom-Singer® CLASSIC™ Indwelling LEIR (Large Esophageal Flange / Increased Resistance)
- ▶ Both voice prosthesis are available as 16 and 20 Fr. version in sizes 6 - 12 IR variant and sizes 4 - 8 LEIR variant
- ▶ With the help of the LARYVOX® INSERTER the voice prosthesis can be placed anterograde quickly and easily:  
 LARYVOX® INSERTER (REF 25900) - for different types of voice prosthesis  
 LARYVOX® INSERTER EXPERT (REF 25910) - for 16 Fr. voice prosthesis  
 LARYVOX® INSERTER EXPERT SMALL (REF 25915) - for 4 mm voice prosthesis

Application with gel capsule technique:

LARYVOX® INSERTER SPECIAL - (REF 25920) - for difficult tracheostoma

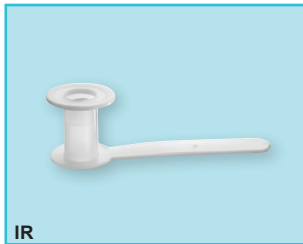

Single patient use  
 Consult instructions for use

### Order information:

<b>Blom-Singer® CLASSIC™ Indwelling IR (Increased Resistance)</b>	<b>16 Fr.</b>	SIZE 6 -12	<b>REF IN16XX-IR</b>
	<b>20 Fr.</b>	SIZE 6 -12	<b>REF IN20XX-IR</b>
<b>Blom-Singer® CLASSIC™ Indwelling LEIR (Large Esophageal Flange / Increased Resistance)</b>	<b>16 Fr.</b>	SIZE 4 - 8	<b>REF IN16XX-LEIR</b>
	<b>20 Fr.</b>	SIZE 4 - 8	<b>REF IN20XX-LEIR</b>
<b>LARYVOX® INSERTER</b>			<b>REF 25900</b>
<b>LARYVOX® INSERTER EXPERT</b>			<b>REF 25910</b>
<b>LARYVOX® INSERTER EXPERT SMALL</b>			<b>REF 25915</b>
<b>LARYVOX® INSERTER SPECIAL</b>			<b>REF 25920</b>