


LARYVOX® TRACHEOSPREIZER

- ▶ The tracheal dilator is used for the dilation of the tracheostoma of tracheostomized and laryngectomized patients
- ▶ To be used for all surgical tracheostomies and as an auxiliary tool for percutaneous techniques
- ▶ Also to be used during cannula changing and for emergency interaction, e.g. decannulation
- ▶ Made of stainless steel
- ▶ The special shape of the rounded tips minimizes the risk of injury during insertion
- ▶ The three spreadable wings allow the opening of the dilator in all directions and thereby ensure an even circular dilation of the tracheostoma
- ▶ Ergonomic design and handling: The dilator opens when the handles are pressed together. This is conform to the natural ergonomics of manual handling and provides a more sophisticated and fine-tuned control of the dilation process.
- ▶ Additional PED variant available (pediatric size)



Consult instructions for use

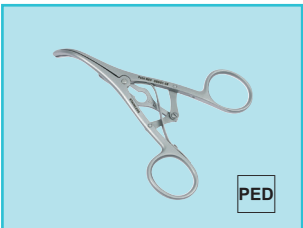

CE marking

Order information:

LARYVOX® Tracheospreizer (adult size)

REF 35500

LARYVOX® Tracheospreizer PED (pediatric size)

REF 35501