
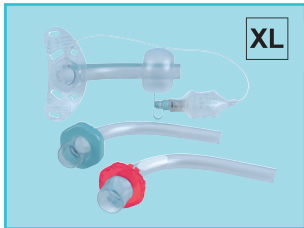


DURACUFF® UNI-VARIO LINGO

1 ICV + 1 ICFU

- ▶ Tracheostomy tube with two inner cannulas (IC), one with a rotating 15 mm swivel connector (ICV / VARIO) to buffer external forces like e.g. tension and torque forces exerted by the ventilation tubing – one sieve-fenestrated with a 15 mm connector (ICFU / UNI) to connect a speaking valve
- ▶ With low-pressure cuff
- ▶ With DURACUFF® Clip for locking the inner cannula into the outer
- ▶ Soft, pliable, lightweight
- ▶ Transparent, inconspicuous design
- ▶ Set includes broad and adjustable tube holder (OPTIFLAUSCH®), obturator and product ID card
- ▶ Three lengths: Standard length, medium (M) and extra-long (XL)

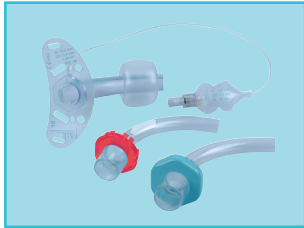

STERILE EO Sterilized using ethylene oxide

No flame Do not resterilize

No sun Keep away from sunlight

Single patient use Single patient use

Read manual Consult instructions for use

Rx only CAUTION: United States Federal law restricts this device to sale, distribution and use by or on order of a physician or a licensed practitioner (in the US only).

CE 0482 CE marking identification number of the notified body

Product made in Germany

Order information:

DURACUFF® VARIO-KOMBI LINGO 1 ICV + 1 ICFK, SIZE 7-12

DURACUFF® VARIO-KOMBI LINGO M 1 ICV + 1 ICFK, SIZE 7-12

DURACUFF® VARIO-KOMBI LINGO XL 1 ICV + 1 ICFK, SIZE 7-12

When ordering, please state the size after the item number!



Size	O.D. neck flange mm	O.D. Tip mm	I.D. mm O-cannula	I.D. mm I-cannula	Length mm standard	Length mm medium (M)	Length mm extra long (XL)
7	10.2	8.8	7.0	5.5	65	75	100
8	11.2	9.8	8.0	6.5	70	80	105
9	12.4	11.2	9.0	7.5	75	85	105
10	13.4	12.2	10.0	8.5	85	95	110
11	14.4	13.2	11.0	9.5	87	97	115
12	15.4	14.2	12.0	10.2	90	100	120

O.D. Tip = Outer diameter of the outer tube at the tip of the cannula · O.D. neck flange = Outer diameter behind the neck flange · I.D. = Inner diameter at the tip of the cannula · Length = Length over the outer curve