


LARYVOX® TOUCH HME HIGHFLOW

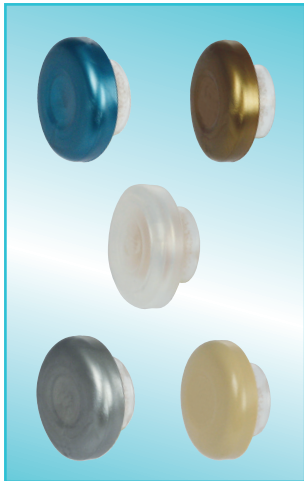

- ▶ Heat and Moisture Exchanger with low-pressure finger-occlusion facility
- ▶ Combines the advantages of an HME with hygienic finger occlusion and an appealing design
- ▶ HighFlow HME with low respiratory resistance and particularly suitable in the case of strong physical exertion and for new HME users
- ▶ The HME is covered by a sealing plate to assure plane occlusion of the HME during speech
- ▶ The sealing plate closes even at a very low pressure of the finger, even if the pressure is not in a perfectly centric position
- ▶ The sealing plate is covered by a large plastic cover to exclude any potential contamination of the finger by secretions or vice versa any potential contamination of the HME by microbial contamination of the finger
- ▶ The inhaled/exhaled air enters/escapes through lateral openings positioned at the side of the HME
- ▶ Therefore no adhesion of bibs or clothing to the HME during inhalation
- ▶ HME with 22 mm adapter (KOMBI) to be used with all LARYVOX® Tape Base Plates or other Base Plate Systems with 22 mm adapter
- ▶ LARYVOX® TOUCH HME HighFlow available in 5 different colours

Order information:

LARYVOX® TOUCH HME HIGHFLOW - blue metallic, Pack of 30

REF 49850-01

LARYVOX® TOUCH HME HIGHFLOW - gold, Pack of 30

REF 49850-02

LARYVOX® TOUCH HME HIGHFLOW - silver, Pack of 30






REF 49850-03

LARYVOX® TOUCH HME HIGHFLOW - beige, Pack of 30

REF 49850-04

LARYVOX® TOUCH HME HIGHFLOW - clear, Pack of 30

REF 49850-13

-  Do not reuse
-  Keep dry
-  Keep away from sunlight
-  Consult instructions for use
-  CE marking