


OPTIBRUSH® SET

Set for Cleaning of Tracheostomy Tubes

- ▶ OPTIBRUSH® CONT Cleaning Box (REF 31200) with a large diameter to allow cleaning of standard tracheostomy tubes. The Cleaning Box contains a removable sieve insert which allows the tracheostomy tube to be gently agitated in the cleaning solution.
- ▶ OPTIBRUSH® CLEAN Cleaning Powder (REF 31110 / 100 g): including measuring spoon for easy and accurate preparation of the cleaning solution; simply dissolve in warm tap water. For a thorough but gentle removal of blood and secretion from tracheostomy tubes.
- ▶ OPTIBRUSH® Cleaning Brushes (REF 31850-08 and REF 31850-10): 2 pieces each of the two most common sizes 08 and 10, i.e. altogether 4 pieces. OPTIBRUSH® Cleaning Brushes are flexible and have a rounded tip to avoid scratching of the soft tracheostomy tube material. The OPTIBRUSH® Cleaning Brushes are single patient use products, i.e. they can be re-used by any individual patient

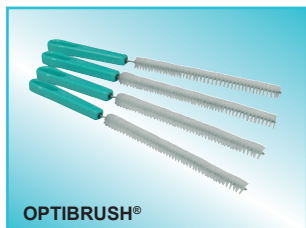

All components of this kit are separately available, cf. detailed individual data sheets

Please refer to the pictograms in the respective instructions of each product.

Order information:

OPTIBRUSH® SET

REF 31011