


## SPIRAFLEX® SHORT UNI LINGO-PHON CUFF ICX, ICFX, SV

- ▶ Short tracheostomy tube with sieve fenestration (LINGO), low pressure cuff (CUFF) and two flexible inner cannulas (ICX, ICFX)
- ▶ Suitable for long-term care
- ▶ Very soft and flexible. The stability is enhanced by the integrated metal spiral.
- ▶ Includes two different flexible inner cannulas one fenestrated (ICFX) and one non-fenestrated cannula (ICX) to facilitate easy tube cleaning and care
- ▶ The flexibility of the outer cannula is not impaired by the use of the inner cannulas
- ▶ Thermoflexible material which softens at body temperature
- ▶ Large, flexible neck flange to reduce skin irritations and pressure sores at the tracheostoma
- ▶ With low pressure cuff (CUFF)
- ▶ Includes two flexible inner cannulas, one fenestrated (ICFX) and one non-fenestrated (ICX), a speaking valve, a cough protection cap, a 15 mm connector (UNI), a closing cap, a tube holder, an obturator and a cannula ID card
- ▶ Both inner cannulas ICX (REF 17516) and ICFX (REF 17517) are also available separately; size to be selected according to the size of the SPIRAFLEX® SHORT outer cannula
- ▶ SPIRAFLEX® SHORT SOLO-ADAPTER PHON (REF 17542) und SPIRAFLEX® SHORT SOLO-ADAPTER UNI (REF 17541) are also available separately

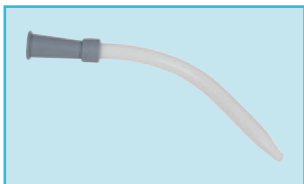

### Order information:

**SPIRAFLEX® SHORT UNI LINGO-PHON CUFF 1 ICX + 1 ICFX + SV, SIZE 6-11**



**REF 17509**

**SPIRAFLEX® SHORT SOLO-ADAPTER UNI**

**REF 17541**

**SPIRAFLEX® SHORT SOLO-ADAPTER PHON**



**REF 17542**

**SPIRAFLEX® SHORT ICX, SIZE 6-11, pack of 3**



**REF 17516**

**SPIRAFLEX® SHORT ICFX, SIZE 6-11, pack of 3**



**REF 17517**

When ordering, please state the size after the item number!

Size	Cuff (Ø) mm	O.D. neck flange mm	O.D. Tip mm	I.D. I mm	I.D. O mm	Length mm	REF
6	20	9.4	9.4	4.5	6.0	61	<b>17509-06</b>
7	22	10.4	10.4	5.5	7.0	67	<b>17509-07</b>
8	24	11.4	11.4	6.5	8.0	77	<b>17509-08</b>
9	25	12.4	12.4	7.5	9.0	77	<b>17509-09</b>
10	26	13.4	13.4	8.5	10.0	88	<b>17509-10</b>
11	28	14.4	14.4	9.5	11.0	88	<b>17509-11</b>



Sterilized using ethylene oxid



Do not resterilize



Keep away from sunlight



Single patient use



Consult instructions for use



CE marking identification number of the notified body

**Product made in Germany**

O.D. Shield = outer diameter behind the neck flange · O.D. Tip = outer diameter at the tube tip · I.D. I = inner diameter at the tip of the inner cannula · I.D. O = inner diameter at the tip of the outer cannula · O-tube = outer cannula · I-tube = inner cannula · Length = length over the outer curve