


# TRACHEOSILC® VARIO

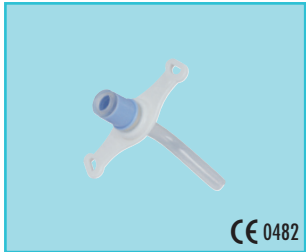

- ▶ Tracheostomy tube without CUFF
- ▶ Also available with CUFF (REF 18830) and CUFF with integrated suction line (SUCTION) to remove secretions above the CUFF (REF 18840 on separate data sheets)
- ▶ Soft, flexible neck-flange with key tube identification data clearly marked
- ▶ Clear and transparent tube shaft for easy handling during cleaning
- ▶ With rotating 15 mm connector (VARIO) to buffer torque and tension caused by tractive forces e.g. from ventilation equipment
- ▶ Rounded cannula tip to prevent irritation of the mucosa in the trachea
- ▶ Set includes tube holder, obturator and product ID card

## Order information:

**TRACHEOSILC® VARIO**

SIZE 5-9.5



**REF 18820**

When ordering, please state the size after the item number!

Size	O.D. neck flange mm	O.D. Tip mm	Length mm	REF
5.0	7.3	7.3	57	<b>18820-05</b>
6.0	8.7	8.7	63	<b>18820-06</b>
7.0	10.0	10.0	71	<b>18820-07</b>
7.5	10.7	10.7	73	<b>18820-075</b>
8.0	11.0	11.0	75	<b>18820-08</b>
8.5	11.7	11.7	78	<b>18820-085</b>
9.0	12.3	12.3	80	<b>18820-09</b>
9.5	13.3	13.3	83	<b>18820-095</b>

O.D. neck flange = outer diameter behind the neck flange · O.D. Tip = outer diameter at the tube tip · Length = length over the outer curve